


the **AQUA-C**
hydration system for **HypoDermoClysis**
...in the 21st Century

"The simple and effective alternative to IV hydration therapy."

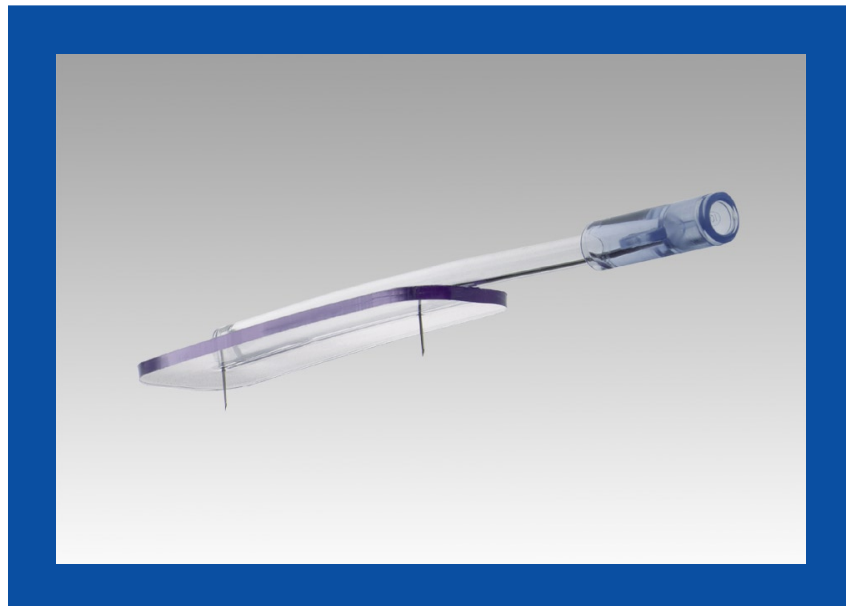

Safe, Simple, Easy to Administer

With the Aqua-C Hydration System, hypodermoclysis is now simple, easy to administer, and relatively painless for the patient.

Hypodermoclysis, the infusion of fluid into the subcutaneous space, is an easy and convenient way of providing hydration for long-term care patients where venous access can be difficult or is unacceptable.

The Aqua-C is easier and less distressing to administer than IV hydration, more suitable for home use, and often eliminates

the need for hospitalization to find a vein. The Aqua-C's smaller 27 gauge x 6mm long needles enhance patient comfort while multiple needles promote even fluid dispersion.

Hypodermoclysis does not cause thrombophlebitis and it is less likely than IV administration to cause pulmonary edema or fluid overload.

The Aqua-C is the perfect solution where IV access is not an option. Call Norfolk Medical for your free sample.

PATENT PENDING

tel: 847-674-7075 • fax: 847-674-7066
info@norfolkmedical.com • www.norfolkmedical.com

Norfolk  Medical

PRODUCT INFORMATION

For your convenience, all of the supplies you need for the safe and simple delivery of fluids to the dehydrated, long-term care patient, are supplied in a sterile pouch.



Aqua-C Hydration System Product Specification Information

- HDC-2500** *High Flow Aqua-C Hydration System* 125cc/hr., Gravity flow, 36"Height
- 2 x (25 gauge x 6mm) needles offset 1 1/2" on a flexible base with adhesive backing, with a female luer, and 1" tubing
 - 82" 20 drop IV administration set with universal spike and flow regulator
 - prep swabstick, alcohol swabs, latex free gloves and transparent dressing
- HDC-2700** *Low Flow Aqua-C Hydration System* 80 cc/hr., Gravity flow, 36"Height
- 2 x (27 gauge x 6mm) needles offset 1 1/2" on a flexible base with adhesive backing, with a female luer, and 1" tubing
 - 82" 20 drop IV administration set with universal spike and flow regulator
 - prep swabstick, alcohol swabs, latex free gloves and transparent dressing

Aqua-C Hydration Needle Set Product Specification Information

- HDC-2521** *High Flow Aqua-C Hydration Needle Set* 125cc/hr., Gravity flow, 36"Height
- 2 x (25 gauge x 6mm) needles offset 1 1/2" on a flexible base with adhesive backing, with a female luer, and 1" tubing
- HDC-2721** *Low Flow Aqua-C Hydration Needle Set* 80 cc/hr., Gravity flow, 36"Height
- 2 x (27 gauge x 6mm) needles offset 1 1/2" on a flexible base with adhesive backing, with a female luer, and 1" tubing

"HypoDermoClysis.....in the 21st Century"